


CERTUS⁺

FLEX LIQUID DISPENSER

User Manual

CERTUS FLEX

CERTUS FLEX

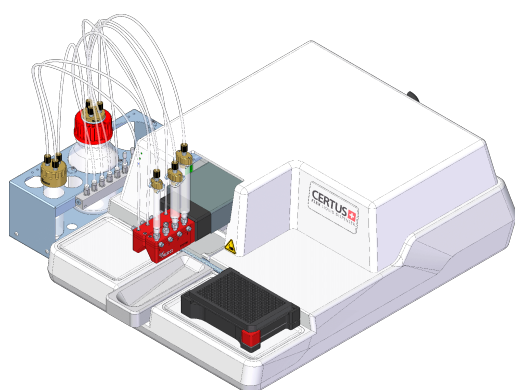

12.0

|

March 2026

GYGER⁺
Fritz Gyger AG

Table of Contents

1. Introduction	0
2. About this Documentation	0
2.1 Additional Documents needed	0
2.2 General	0
2.2.1 Language version	0
2.2.2 Electronic documentation	0
2.3 Copyright	0
2.4 Responsibility for the documentation	0
2.5 Maintenance and storage of documentation	0
2.5.1 Maintenance	0
2.5.2 Quality	0
2.6 Presentation of Information	0
2.6.1 Safety Instructions	0
2.7 Document Layout	0
2.8 Signal Words	0
2.9 Standards and Directives applied	0
2.9.1 2014/35/EU	0
2.9.2 2014/30/EU	0
2.9.3 2006/42/EC	0
2.10 Text Conventions	0
2.10.1 Paths in software	0
2.10.2 Real actions	0
2.11 Presentation example	0
2.11.1 English terms	0
3. Safety	0
3.1 Intended Use	0
3.2 Non-Intended Use	0
3.3 Personnel Qualifications	0
3.3.1 Requirements	0
3.4 Basic Safety Regulations	0
3.4.1 Regulations and Work Safety	0
3.4.2 General Inspections and Maintenance	0
3.4.3 Alterations	0
3.4.4 Power Supply	0

3.4.5 Spare Parts	0
3.4.6 Fluid Media	0
3.4.7 Personal Protective Equipment	0
3.4.8 Dealing with Malfunctions	0
3.4.9 Cleaning, Disinfection, and Sterilization	0
3.4.10 Disposal	0
3.5 Responsibilities	0
3.5.1 Responsibilities of the End User	0
3.5.2 Responsibilities of Operating Personnel	0
3.5.3 Responsibilities of Maintenance Personnel	0
3.6 Safety Instructions	0
3.6.1 General Safety Instructions	0
3.6.2 Statutory Safety Instructions	0
3.6.3 Consequences if Safety Instructions are Not Followed	0
3.7 Special Hazards	0
3.7.1 Electrical Equipment	0
3.7.2 Danger from Moving Parts	0
3.7.3 Danger from Toxic or Highly Reactive Substances	0
3.7.4 Danger from Substances Harmful to Human Health	0
3.8 First Aid in Case of Accident	0
3.8.1 Training Areas	0
4. Transportation and Storage	0
4.1 Transportation – Safety Instructions	0
4.1.1 Risk of crushing	0
4.1.2 Damage due to inadequate transportation	0
4.2 Transportation – Inspection	0
4.3 Packaging	0
4.3.1 Special labelling of the packaging	0
4.3.2 Disposal of packaging materials	0
4.4 Storage	0
5. Technical Specifications - Dimensions	0
5.1 Base Unit	0
5.2 System Integration	0
5.2.1 Mounting Point Specifications	0
5.3 Operating Range	0
5.3.1 Operating Range Dimensions	0

6. Technical Specifications - Device Data	0
6.1 Base Unit – General Data	0
6.1.1 Physical Specifications	0
6.1.2 Axis Performance	0
6.1.3 Electrical Specifications	0
6.2 Ambient Conditions	0
6.3 Compressed Air	0
6.4 Dispensing Head	0
6.5 Fluid and Supply	0
6.6 Well Plates	0
6.7 Dispensing Times (Examples)	0
6.7.1 Performance Data	0
6.8 Micro valve and Dispensing	0
6.9 Regulatory	0
7. Start Up	0
7.1 General Safety Instructions	0
7.2 Preparation	0
7.3 Unpacking	0
7.4 Warranty	0
7.5 CERTUS FLEX – Start-Up	0
7.5.1 Inserting the Waste Tray	0
7.5.2 Connecting the Unit	0
7.5.3 Putting on the Dispensing Head	0
7.5.4 Installation of the Micro Valves	0
7.5.5 Preparing the Bottle Kit	0
7.5.6 Step 1: Cutting tubes to length	0
7.5.7 Step 2: Placing the clamping sleeve	0
7.5.8 Step 3: Inserting the inlet adapter	0
7.5.9 Step 4: Driving tubes into bottle caps	0
7.5.10 Step 5: Installing a suction filter (option)	0
7.5.11 Step 6: Capping the lab bottle	0
7.5.12 Preparing the Syringe Kit	0
7.5.13 Step 1: Cutting the compressed air tube to length	0
7.5.14 Step 2: Placing the clamping sleeve	0
7.5.15 Step 3: Connecting the tube to the adapter	0
7.5.16 Step 4: Screwing the inlet adapter into the syringe	0
7.5.17 Step 5: Installing the adapter head onto the syringe	0

7.5.18	Installing the Bottle/Syringe Kit on the Dispensing Head	0
7.5.19	Refilling the Fluid Container	0
7.6	Refilling bottle kits	0
7.7	Refilling syringe kits	0
8.	Assembly and Operation	0
8.1	Directions of Movement	0
8.2	CERTUS FLEX – Assembly	0
8.3	CERTUS FLEX – Operating Principle	0
8.4	Base Unit	0
8.4.1	Rating Plate and Interfaces	0
8.4.2	Status Indicator	0
8.5	Bottle and Syringe Kit	0
8.5.1	Bottle Kit – Assembly	0
8.5.2	Bottle Kit – Operational Chart	0
8.5.3	Bottle Kit – Bottle Caps	0
8.5.4	Syringe Kit – Assembly	0
8.5.5	Syringe Kit – Operational Chart	0
8.6	Bottle Holder	0
8.7	Dispensing Head	0
8.7.1	Dispensing head – 22°	0
8.7.2	Magnetic stirrer dispensing head	0
8.8	SMLD 300GC micro valve	0
8.8.1	Overview – Dispensing Range	0
8.9	Waste Tray	0
8.10	Well Plate Holder	0
9.	Operation	0
9.1	Safety Instructions for Operators	0
9.2	Actions Prior to Use	0
9.3	Power-Up	0
9.4	Overview - Experiment	0
9.5	Gathering Data	0
9.5.1	Fluid	0
9.5.2	Micro valves	0
9.5.3	Well Plates	0
9.5.4	Standard 96 Well Plate	0
9.5.5	Standard 384 Well Plate	0
9.5.6	Standard 1536 Well Plate	0

9.5.7	Calibrations	0
9.6	Opening an Experiment/Protocol	0
9.7	Creating a New Protocol	0
9.8	Configuring a Dispensing Head	0
9.8.1	Using a Dispensing Head Template	0
9.8.2	Copy the Configuration Template	0
9.8.3	Saving an Active Configuration	0
9.8.4	Defining a Supply Kit	0
9.8.5	Configuration transfer to another dispensing head type	0
9.9	Maintenance	0
9.9.1	Preparation	0
9.10	Calibration	0
9.10.1	Determining Calibration Points	0
9.10.2	Adding customized Calibration Points	0
9.10.3	Completing Metadata	0
9.10.4	Verification using the Graph	0
9.10.5	Volume-Graph	0
9.10.6	Verifying a Calibration	0
9.10.7	Verifying a User-Defined Volume	0
9.11	Defining a Protocol	0
9.11.1	Creating a New Protocol	0
9.11.2	Refilling the Evaporation Curtain	0
9.11.3	Regions	0
9.11.4	Layer	0
9.12	Import CSV file	0
9.12.1	Import CSV-volumes into a new layer	0
9.12.2	Edit CSV volumes	0
9.12.3	CSV format	0
9.13	Conducting an experiment	0
9.14	Interval Purge Loop (IPL)	0
10.	CERTUS CONTROL	0
10.1	Features	0
10.2	Minimum System Requirements	0
10.3	Installation	0
10.3.1	Installation Files	0
10.3.2	Configuring the IP Address	0
10.3.3	Installing the INCO Server	0

10.3.4	Installing CERTUS FLEX Control	0
10.3.5	Stand-Alone Application	0
10.3.6	Integrated Application	0
10.3.7	Additional information	0
10.4	Logging on	0
10.4.1	Logging on as an Administrator	0
10.4.2	Logging onto the Simulation Mode	0
10.5	Overview -- Start Menu	0
10.6	Library	0
10.6.1	Menu bar	0
10.6.2	Fluid	0
10.6.3	Valve	0
10.6.4	Well Plate	0
10.6.5	Evaporation Curtain	0
10.6.6	Calibration	0
10.7	Settings	0
10.7.1	Default Values	0
10.7.2	Protocol Editor	0
10.7.3	Maintenance	0
10.7.4	Device	0
10.8	Help	0
10.9	Users	0
10.9.1	Symbols in use	0
10.9.2	Context menu	0
10.9.3	Adding Users	0
10.10	Protocol Editor	0
10.10.1	Symbols in use	0
10.10.2	Operating Principle	0
10.10.3	Experiment Menu	0
10.10.4	Configuration	0
10.10.5	Template Manager	0
10.11	CERTUS Calibrator	0
10.11.1	Symbols in use	0
10.11.2	Context menu	0
10.11.3	Context menu	0
10.11.4	Configuration Manager	0
10.11.5	Symbols in use	0

10.11.6	Context menu	0
10.11.7	Context menu	0
10.11.8	Calibration Manager	0
10.11.9	Maintenance	0
10.11.10	Live Log	0
10.11.11	Summary	0
10.11.12	Interval Purge Loop (IPL)	0
10.12	Protocol	0
10.12.1	Symbols in use	0
10.12.2	Context menu	0
10.13	Region Editor	0
10.13.1	Add dispense step	0
10.13.2	Synchronous dispense	0
10.13.3	Axial movement	0
10.13.4	Sequence Replicate	0
10.14	Dispensing Functions	0
10.14.1	Constant	0
10.14.2	Linear	0
10.14.3	Factor	0
10.14.4	Discrete	0
10.14.5	Backfill	0
11.	Maintenance	0
11.1	General Safety Instructions	0
11.2	Cleaning, Disinfection, and Sterilization	0
11.2.1	Manufacturer's Cleaning Recommendations	0
11.3	Cleaning Table	0
11.4	Maintenance Works	0
11.4.1	Replacing Spare Parts	0
11.4.2	General Cleaning	0
11.4.3	Cleaning the Micro Valve	0
11.4.4	Removing/Replacing the micro valve	0
11.4.5	Fluid Change – Syringe Kit	0
11.4.6	Fluid Change -Bottle Kit	0
11.4.7	Replacing the Fuse on the Base Unit	0
11.4.8	Removing the Dispensing Head	0
12.	Troubleshooting	0
12.1	General Safety Instructions	0

12.2	What to do in case of a Malfunction	0
12.3	CERTUS FLEX - Hardware	0
12.3.1	CERTUS CONTROL - Protocol Validator	0
12.3.2	CERTUS CONTROL - Maintenance	0
13.	Decommissioning and Disposal	0
13.1	General Safety Instructions	0
13.2	Decommissioning the Unit	0
13.3	Storing the Unit	0
13.4	Disposing Unit	0
14.	References	0
15.	Change Index	0
16.	Decontamination Declaration	0
16.1	Decontamination Declaration for CERTUS FLEX and Micro Valves	0
16.1.1	Customer details	0
16.1.2	Decontamination / Disinfection	0
16.1.3	Declaration	0
17.	Document List	0
18.	WEB CONTROL	0
18.1	WEB CONTROL - Installation	0
18.2	WEB CONTROL - Operation	0
18.2.1	Open the control page in the web browser	0
18.3	WEB CONTROL - Operation	0
18.3.1	Configuration of IP address and port	0
18.3.2	Predefinition of prime / purge volumes	0
18.3.3	Execute commands in the web browser	0
18.4	WEB CONTROL - Error Messages	0